

Plant Hormone abscisic acid ELISA Kit (ABA)

Catalog NO.: RK04822

version: 2.0

This package insert must be read in its entirety
before using this product

Introduction

This ELISA kit used for quantitative determination of ABA in Plant tissue and Other biological fluids. For research use only, and it is highly recommended to read thoroughly of this manual before using the product.

Principle of the Assay

This assay is based on the competitive binding technique in which ABA present in a sample competes with a fixed amount of horseradish peroxidase (HRP)-labeled ABA for sites on a Plant monoclonal antibody. Following a wash to remove excess conjugate and unbound sample, a chromogenic substrate solution is added to the wells to determine the bound enzyme activity. The color development is stopped by acid, and the absorbance is read at 450 nm. The intensity of the color is inversely proportional to the concentration of Cortisol in the sample.

Material Provided & Storage Conditions

Unopened kits can be stored at 2-8°C for 1 year, and opened products must be used within 1 month.

Part	Size	Cat. No.	Storage of opened/reconstituted material
Antibody Coated Plate	8x12	RM97000	Put the unused slats back in the aluminum foil bag with the desiccant and reseal them. They can be stored at 2-8°C for 1 month.
Standard Lyophilized	2 vials	RM97001	It is not recommended to use again after redissolving.
Concentrated HRP- Conjugate Antigen(100x)	1 ×60μL	RM97003	Store at 2-8°C for 1 month *
Standard/Sample Diluent (R1)	1 ×20mL	RM00023	Store at 2-8°C for 1 month *
HRP- Conjugate Antigen Diluent(R2)	1 ×6mL	RM00024	

Wash Buffer(25x)	1 ×20mL	RM00026	Store at 2-8°C for 1 month *
TMB Substrate	1 ×12mL	RM00027	
Stop Solution	1 ×6mL	RM00028	
Plate Sealers	4 Strips		
Specification	1		

Other Supplies Required

1. Microplate reader capable of measuring absorbance at 450 nm, with the correction wavelength set at 630 nm or 570 nm.
2. Pipettes and pipette tips.
3. Deionized or distilled water.
4. Squirt bottle, manifold dispenser, or automated microplate washer.
5. Incubator.
6. Test tubes for dilution of standards and samples.

Precautions

*** FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.**

1. Any variation in diluent, operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.
2. Variations in sample collection, processing, and storage may cause sample value differences.
3. Reagents may be harmful, if ingested, rinse it with an excess amount of tap water.
4. Stop Solution contains strong acid. Wear eye, hand, and face protection.
5. Please perform simple centrifugation to collect the liquid before use.
6. Do not mix or substitute reagents with those from other lots or other sources.
7. Adequate mixing is particularly important for good result. Use a mini-vortexer at the lowest frequency.
8. Mix the sample and all components in the kits adequately and use clean plastic container to prepare all diluents.
9. Both the sample and standard should be assayed in duplicate, and reagents should be added in sequence in

- accordance with the requirement of the specification.
10. Reuse of dissolved standard is not recommended.
 11. The kit should not be used beyond the expiration date on the kit label.
 12. The kit should be away from light when it is stored or incubated.
 13. To reduce the likelihood of blood-borne transmission of infectious agents, handle all serum, plasma, and other biological fluids in accordance with NCCLS regulations.
 14. To avoid cross contamination, please use disposable pipette tips.
 15. Please prepare all the kit components according to the Specification. If the kits will be used several times, please seal the rest strips and preserve with desiccants. Do use up within 2 months.
 16. This assay is designed to eliminate interference by other factors present in biological samples.
 17. Until all factors have been tested in this assay, the possibility of interference cannot be excluded.
 18. The 48T kit is also suitable for the specification.

Sample Collection & Storage

The sample collection and storage conditions listed below are intended as general guidelines. Sample stability has not been evaluated.

Plant tissues (1): Take 500 mg clean and air-dry plant tissue. Put them in the mortar, mix with liquid nitrogen and grind quickly until it becomes fine powders, add 500 μL 1X PBS, mix. Then after 10000rpm centrifugation for 5 min, take 50 μL of the supernatant for analysis. According to the result of pretest, determine whether it is needed for dilution with the 1X PBS

Plant tissues (2): Clean and grind plant tissues. Add 500 mg plant tissues to 500 μL 1X PBS. Put them on ice for a half-hour. Then homogenized them with a homogenizer, or with a ultrasonic disrupter until there is no significant tissue mass. Then after 10000rpm centrifugation for 5 min, take 50 μL of the supernatant for analysis. According to the result of pretest, determine whether it is needed for dilution with the 1X PBS..

Other biological fluids : Centrifuge samples for 20 minutes at $1,000 \times g$. Collect the supernatants and assay immediately or store samples in aliquot at -20°C or -80°C for later use. Avoid

repeated freeze-thaw cycles.

Note: Samples to be used within 5 days may be stored at 2-8 °C, other wise samples must be stored at -20°C (\leq 1 month) or -80°C (\leq 2 months) to avoid loss of bioactivity and contamination.

SAMPLE PREPARATION

Plant tissues samples require no dilution before test. The recommended dilution factor is for reference only. The optimal dilution factor should be determined by users according to their particular experiments.


Reagent Preparation

Bring all reagents to room temperature before use. If crystals have formed in the concentrate, Bring the reagent to room temperature and mix gently until the crystals have completely dissolved.

Standard - Reconstitute the Standard Lyophilized with 1.0 mL

Standard/Sample Diluent(R1). This reconstitution produces a stock solution of 100ng/mL. Mix the standard to ensure complete reconstitution and allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions.

Use the 100ng/mL standard stock to produce a dilution series (below) with Standard/Sample Diluent(R1). Mix each tube thoroughly and change pipette tips between each transfer (recommended concentration for standard curve: 100, 50, 25, 12.5, 6.25, 3.12, 1.56, 0ng/mL). Use diluted standards within 60 minutes of preparation.

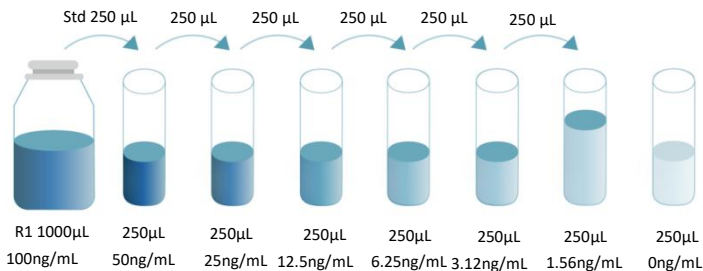

Working HPR-Conjugate Antigen - Dilute 1:100 of Concentrated HRP-Conjugate Antigen (100x) with HRP-Conjugate Antigen Diluent (R2) before use. For example: Add 10 µ L of

Concentrated HRP-Conjugate Antigen (100x) to 990 μ L
HRP-Conjugate Antigen Diluent (R2) to prepare 1000 μ L
Working HRP-Conjugate Antigen Buffer.

Wash Buffer - If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved. Dilute 1:25 with double distilled or deionized water before use, for example : Add 16 mL of Wash Buffer Concentrate to 384 mL of deionized or distilled water to prepare 400 mL of Wash Buffer.

Assay Procedure

Bring all reagents and samples to room temperature before use. It is recommended that all standards, controls, and samples be assayed in duplicate.

1. Prepare all reagents, working standards, and samples as directed in the previous sections.
2. Remove excess microplate strips from the plate frame, return them to the foil pouch containing the desiccant pack, and reseal.
3. Add 50 μ L Standard/sample Diluent (R1) in a blank well.

4. Add 50 μ L different concentration of standard or sample in other wells.. Record the plate layout of standards and sample assay.
5. Add 50 μ L Working HRP-Conjugate Antigen in each wells (**Note: Do not add to Blank!**) , and cover with new adhesive sealer provided.Cover with the adhesive strip provided. Incubate for 1 hours at 37°C.
6. During the incubation, turn on the microplate reader to warm up for 30 minutes before measuring.
7. Add wash buffer 250 μ L/well, aspirate each well after holding 1-2 minutes, repeating the process four times for a total of five washes.
8. Add 100 μ L TMB Substrate to each well. Incubate for 15-20 minutes at 37°C. Protect from light.
9. Add Stop Solution (50 μ L/well), determine the optical density of each well within 5 minutes, using a Microplate reader set to 450 nm. If wavelength correction is available, set to 570 nm or 630 nm. If wavelength correction is not available, subtract readings at 570 nm or 630 nm from the readings at 450 nm. This subtraction will correct for optical imperfections in the plate. Readings

made directly at 450 nm without correction may cause higher value and less accurate result.

Assay Procedure Summary

Prepare the standard and reagents



Add 50 μ L of standards or test samples to each well



Add 50 μ L Working HRP-Conjugate Antigen (**Note: Do not add to Blank!**)

Incubate for 1hours at 37°C, then wash 5 times



Add 100 μ L Substrate Solution

Incubate for 15-20 min at 37°C under dark condition



Add 50 μ L Stop Solution




Detect the optical density within 5 minutes under 450nm.

Correction Wavelength set at 570nm or 630nm

Calculation of Results

1. Average the duplicate readings for each standard, control and sample, and subtract the average zero standard optical density (O.D.).
2. Create a standard curve by reducing the data using computer software capable of generating a four-parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the Y-axis against the concentration on the X-axis and draw a best fit curve through the points on a log/log graph. The data may be linearized by plotting the log of the ABA concentrations versus the log of the O.D. on a linear scale, and the best fit line can be determined by regression analysis.
3. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

Typical Data

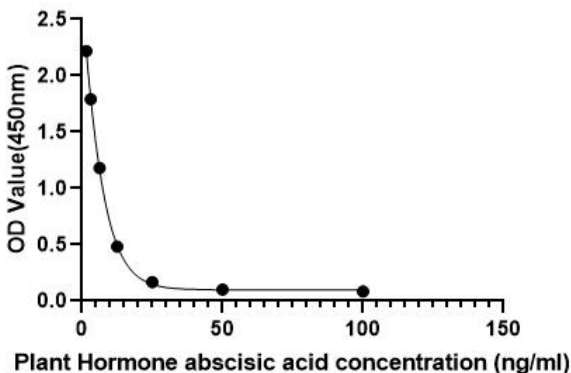

The standard curves are provided for demonstration only. A standard curve should be generated for each set of ABA assayed.

Detection Range

1.56-100ng/mL

Sensitivity

The minimum detectable dose (MDD) of ABA typically less than

0.78ng/mL. The MDD was determined by adding two standard deviations to the mean optical density value of twenty zero standard replicates and calculating the corresponding concentration.

Specificity

This assay has high sensitivity and excellent specificity for detection of ABA .

No significant cross-reactivity or interference between ABA and analogues was observed.

Note:

Limited by current skills and knowledge, it is impossible for us to complete the cross-reactivity detection between ABA and all the analogues, therefore, cross reaction may still exist.

Precision

Intra-plate Precision

3 samples with low, middle and high level ABA were tested 20 times on one plate, respectively.

Intra-Assay: CV<10%

Inter-plate Precision

3 samples with low, middle and high level ABA were tested on 3 different plates, 20 replicates in each plate.

Inter-Assay: CV<15%

Recovery

Matrices listed below were spiked with certain level of ABA and the recovery Percentages were calculated by comparing the measured value to the expected amount of ABA in samples.

Sample	Average Recovery (%)	Range (%)
Cell Culture Media(n=5)	88	83-103
Serum(n=5)	84	81-90

Linearity

The linearity of the kit was assayed by testing samples spiked with appropriate concentration of ABA and their serial dilutions. The results were demonstrated by the percentage of calculated concentration to the expected.

/	/	Cell Culture Media(n=5)	Serum(n=5)
1:2	Average of Expected (%)	92	90
	Range (%)	85-102	80-94
1:4	Average of Expected (%)	100	86
	Range (%)	92-113	87-97
1:8	Average of Expected (%)	96	85
	Range (%)	95-100	80-97
1:16	Average of Expected (%)	103	97
	Range (%)	90-110	84-106

Trouble Shooting

Problem	Possible Cause	Solution
High Background	Insufficient washing	Sufficiently wash plates as required. Ensure appropriate duration and number of washes. Ensure appropriate volume of wash buffer in each well.
	Incorrect incubation procedure	Check whether the duration and temperature of incubation are set up as required.
	Cross-contamination of samples and reagents	Be careful of the operations that could cause cross-contamination. Use fresh reagents and repeat the tests.
No signal or weak signal	Incorrect use of reagents	Check the concentration and dilution ratio of reagents. Make sure to use reagents in proper order.
	Incorrect use of microplate reader	Warm the reader up before use. Make sure to set up appropriate main wavelength and correction wavelength.
	Insufficient colour reaction time	Optimum duration of colour reaction should be limited to 15-25 minutes.

	Read too late after stopping the colour reaction	Read the plate in 5 minutes after stopping the reaction.
	Matrix effect of samples	Use positive control.
Too much signal	Contamination of TMB substrate	Check if TMB substrate solution turns blue. Use new TMB substrate solution.
	Plate sealers reused	Use a fresh new sealer in each step of experiments.
	Protein concentration in sample is too high	Do pre-test and dilute samples in optimum dilution ratio.
Poor Duplicates	Uneven addition of samples	Check the pipette. Periodically calibrate the pipette.
	Impurities and precipitates in samples	Centrifuge samples before use.
	Inadequate mixing of reagents	Mix all samples and reagents well before loading.

*For research purposes only. Not for therapeutic or diagnostic purposes.