

IP6K1 Rabbit pAb

Catalog NO.: A12216

Basic Information

Observed MW

57kDa

Calculated MW

50kDa

Category

Polyclonal Antibodies

Applications

WB, IHC-P, ELISA

Cross-Reactivity

Human, Mouse, Rat

Conjugation

Unconjugated

Background

This gene encodes a member of the inositol phosphokinase family. The encoded protein may be responsible for the conversion of inositol hexakisphosphate (InsP6) to diphosphoinositol pentakisphosphate (InsP7/PP-InsP5). It may also convert 1,3,4,5,6-pentakisphosphate (InsP5) to PP-InsP4. Alternatively spliced transcript variants have been described.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:50 - 1:200

ELISA

Recommended starting concentration is 1 µg/mL.
Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID

9807

Swiss Prot

Q92551

Immunogen

Recombinant protein (or fragment).This information is considered to be commercially sensitive.

Synonyms

PiUS; IHPK1; IP6K1


Contact

Product Information


☒		order@abclonal.com
☒		service@abclonal.com
☒		www.abclonal.com

Source	Isotype	Purification
Rabbit	IgG	Affinity purification
Storage		
Store at -20°C. Avoid freeze / thaw cycles.		
Buffer: PBS with 0.01% thimerosal,50% glycerol,pH7.3.		


Validation Data

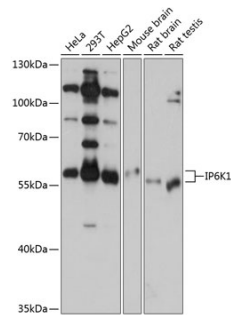

Western blot analysis of various lysates using IP6K1 Rabbit pAb (A12216) at 1:1000 dilution.
 Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.
 Lysates/proteins: 25µg per lane.
 Blocking buffer: 3% nonfat dry milk in TBST.
 Detection: ECL Basic Kit (RM00020).
 Exposure time: 10s.

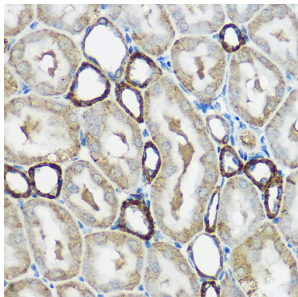

Immunohistochemistry analysis of paraffin-embedded Rat kidney using IP6K1 Rabbit pAb (A12216) at dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

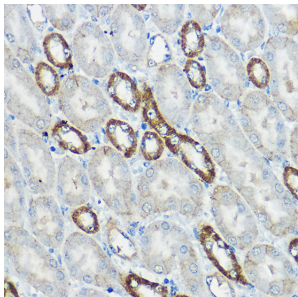

Immunohistochemistry analysis of paraffin-embedded Mouse kidney using IP6K1 Rabbit pAb (A12216) at dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate buffer (pH 6.0) prior to IHC staining.

