

pan-TriMethyl-lysine Rabbit pAb

Catalog NO.: A18292 **1 Publications**

Basic Information

Observed MW

18-55kDa

Calculated MW

Category

Methylated Antibodies

Applications

WB, ELISA

Cross-Reactivity

Human, Mouse, Rat

Background

Methylation of lysine residues is a common regulatory post-translational modification (PTM) that results in the mono-, di-, or tri-methylation of lysine at ε-amine groups by protein lysine methyltransferases (PKMTs).The post-translational ε-amino lysine methylated proteins is an important reversible modification which plays a vital role in the regulation of many cellular processes including chromatin dynamics and gene transcription. Methylation of lysine residues is modulated by specific counteractive enzymes including lysine methylases (KMTs) and demethylases (KDMs). Lysine trimethylation occurs in both histones and non-histone substrates. It has become promising targets for discovery of anti-cancer drugs.

Recommended Dilutions

WB 1:500 - 1:2000

ELISA

Recommended starting concentration is 1 µg/mL.
Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID **Swiss Prot**

Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

Synonyms


Contact

Product Information

☒		order@abclonal.com
☒		service@abclonal.com
☒		www.abclonal.com

Source	Isotype	Purification
Rabbit	IgG	Affinity purification
Storage		
Store at -20°C. Avoid freeze / thaw cycles.		
Buffer: PBS containing 50% glycerol, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.		

Validation Data

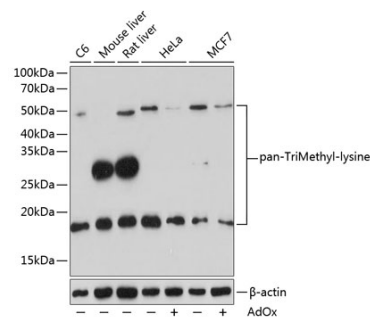

Western blot analysis of various lysates using pan-TriMethyl-lysine pAb (A18292) at 1:500 dilution. HeLa and MCF7 cells were treated by ADOX (100 μM) for 24 hours .

Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25μg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Enhanced Kit (RM00021).

Exposure time: 3min.

