

E-Cadherin Rabbit mAb

Catalog NO.: A20798 **Recombinant** **115 Publications**

Basic Information

Observed MW

135 kDa

Calculated MW

97 kDa

Category

Monoclonal Antibodies

Applications

WB, IF/ICC, IF-F, IF-P, IHC-P, mIHC, ELISA

Cross-Reactivity

Human, Mouse, Rat

Conjugation

Unconjugated

Clone number

ARC51012

Background

This gene encodes a classical cadherin of the cadherin superfamily. Alternative splicing results in multiple transcript variants, at least one of which encodes a preproprotein that is proteolytically processed to generate the mature glycoprotein. This calcium-dependent cell-cell adhesion protein is comprised of five extracellular cadherin repeats, a transmembrane region and a highly conserved cytoplasmic tail. Mutations in this gene are correlated with gastric, breast, colorectal, thyroid and ovarian cancer. Loss of function of this gene is thought to contribute to cancer progression by increasing proliferation, invasion, and/or metastasis. The ectodomain of this protein mediates bacterial adhesion to mammalian cells and the cytoplasmic domain is required for internalization. This gene is present in a gene cluster with other members of the cadherin family on chromosome 16.

Recommended Dilutions

WB	1:1000 - 1:8000
IF/ICC	1:100 - 1:1000
IF-F	1:200 - 1:600
IF-P	1:100 - 1:1000
IHC-P	1:400 - 1:2000
mIHC	1:200 - 1:800

ELISA

Recommended starting concentration is 1 µg/mL.
Please optimize the concentration based on your specific assay requirements.

Contact

Immunogen Information

Gene ID	Swiss Prot
999	P12830

Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

Synonyms

UVO; CDHE; ECAD; LCAM; Arc-1; BCDS1; CD324; E-Cadherin

Product Information

☒		order@abclonal.com
☒		service@abclonal.com
☒		www.abclonal.com

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol and 0.05% BSA, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

Validation Data

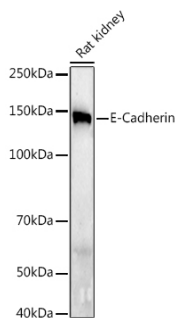

Western blot analysis of lysates from Rat kidney using E-Cadherin Rabbit mAb (A20798) at 1:1000 dilution incubated overnight at 4°C.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25 µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Enhanced Kit (RM00021).

Exposure time: 180s.

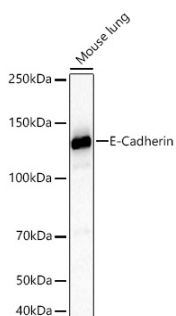

Western blot analysis of lysates from Mouse lung using E-Cadherin Rabbit mAb (A20798) at 1:2000 dilution incubated overnight at 4°C.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25 µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 45s.

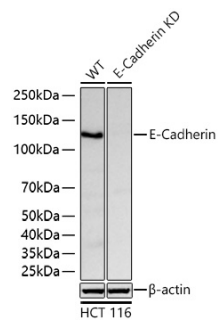

Western blot analysis of lysates from wild type (WT) and E-Cadherin knockdown (KD) HCT 116 cells using E-Cadherin

Rabbit mAb (A20798) at 1:4000 dilution incubated overnight at 4°C.

Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25 µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 30 s.

Validation Data

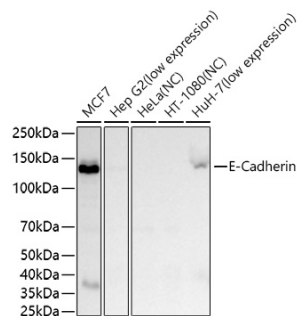

Western blot analysis of various lysates using E-Cadherin Rabbit mAb (A20798) at 1:4000 dilution incubated overnight at 4°C.

Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25 µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.


Detection: ECL Basic Kit (RM00020).

Negative control (NC): HeLa, HT-1080.


Exposure time: 30 s.

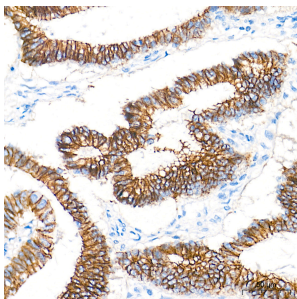


Immunohistochemistry analysis of paraffin-embedded Human colon carcinoma tissue using E-Cadherin Rabbit mAb (A20798) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

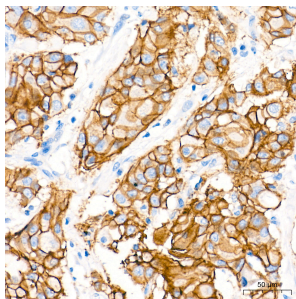


Immunohistochemistry analysis of paraffin-embedded Human breast cancer tissue using E-Cadherin Rabbit mAb (A20798) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

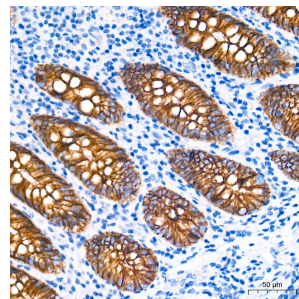


Immunohistochemistry analysis of paraffin-embedded Human colon tissue using E-Cadherin Rabbit mAb (A20798) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

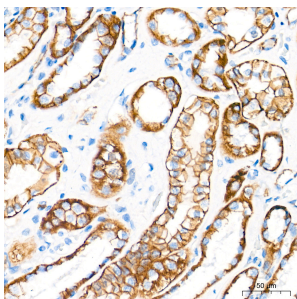

Immunohistochemistry analysis of paraffin-embedded Human kidney tissue using E-Cadherin Rabbit mAb (A20798) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M

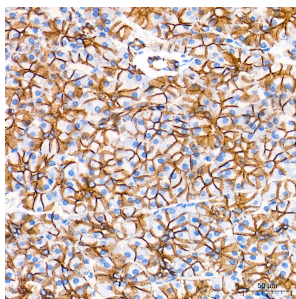

Immunohistochemistry analysis of paraffin-embedded Human pancreas tissue using E-Cadherin Rabbit mAb (A20798) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M

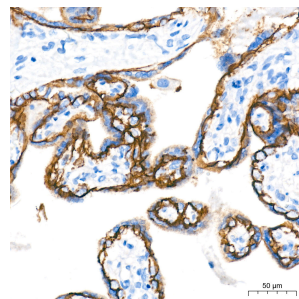


Immunohistochemistry analysis of paraffin-embedded Human placenta tissue using E-Cadherin Rabbit mAb (A20798) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M

Tris-EDTA Buffer(pH 9.0) prior to IHC staining.


Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

Tris-EDTA Buffer(pH 9.0) prior to IHC staining.


Validation Data

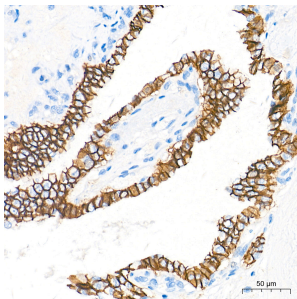


Immunohistochemistry analysis of paraffin-embedded Human thyroid cancer tissue using E-Cadherin Rabbit mAb (A20798) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

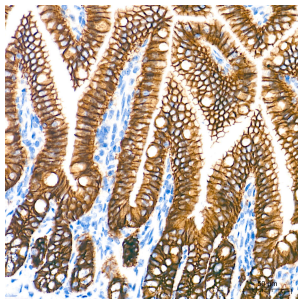


Immunohistochemistry analysis of paraffin-embedded Mouse intestin tissue using E-Cadherin Rabbit mAb (A20798) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

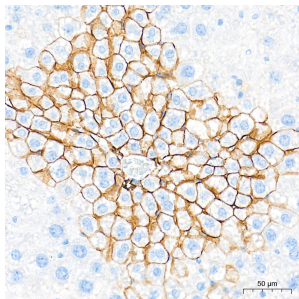


Immunohistochemistry analysis of paraffin-embedded Mouse liver tissue using E-Cadherin Rabbit mAb (A20798) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

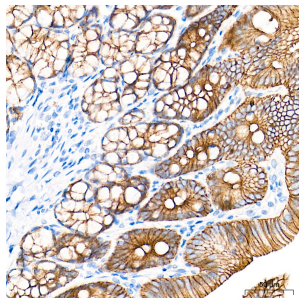


Immunohistochemistry analysis of paraffin-embedded Rat colon tissue using E-Cadherin Rabbit mAb (A20798) at a dilution of 1:500 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer(pH 9.0) prior to IHC staining.

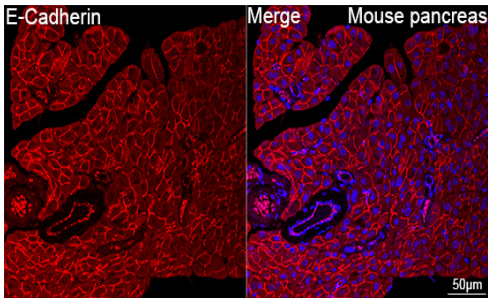


Confocal imaging of paraffin-embedded Mouse pancreas tissue using E-Cadherin Rabbit mAb (A20798, dilution 1:100) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

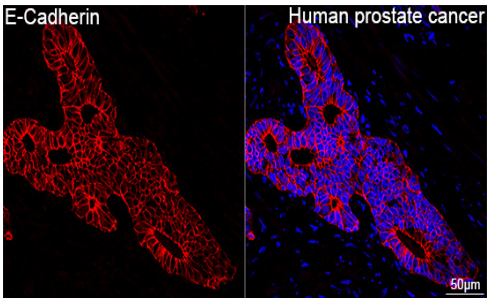


Confocal imaging of paraffin-embedded Human prostate cancer tissue using E-Cadherin Rabbit mAb (A20798, dilution 1:100) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

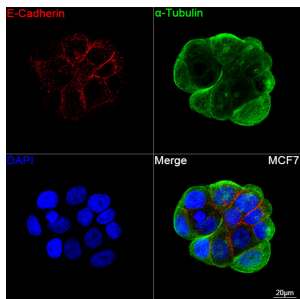

Confocal imaging of MCF7 cells using E-Cadherin Rabbit mAb (A20798, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were

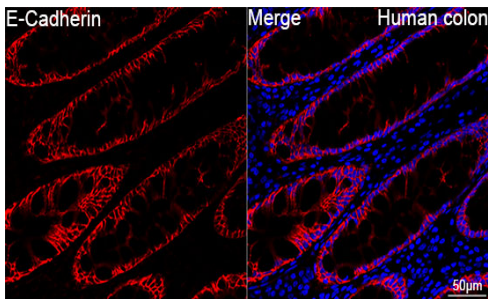

Confocal imaging of paraffin-embedded Humanc colon tissue using E-Cadherin Rabbit mAb (A20798, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI

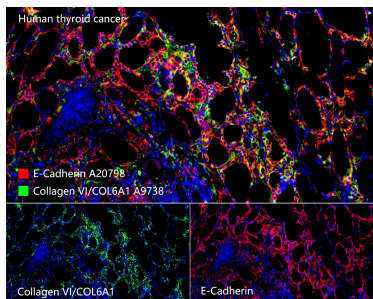

The multiplex IHC analysis on paraffin-embedded Human thyroid cancer tissue using the following specific primary antibodies and tyramide signal amplification (TSA) reagents (RK05903) : Collagen

counterstained with α -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

Antibody | Protein | ELISA Kit | Enzyme | HGS | Service


was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

VI/COL6A1 Rabbit mAb (A9738, 1:100) with TSA-TYR-520 (Green), and E-Cadherin Rabbit mAb (A20798, 1:500) with TSA-TYR-570 (Red). DAPI (Blue) was used for nuclear staining. Prior to multiplex IHC staining, high-pressure antigen retrieval was performed using 0.01M citrate buffer at pH 6.0. The analysis was completed using a 40x objective lens.


For research use only. Not for therapeutic or diagnostic purposes.

Please visit <https://abclonal.com> for a complete list of recommended products.


Validation Data

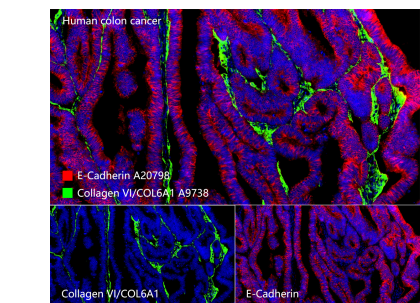


The multiplex IHC analysis on paraffin-embedded Human colon cancer tissue using the following specific primary antibodies and tyramide signal amplification (TSA) reagents (RK05903) : Collagen VI/COL6A1 Rabbit mAb (A9738, 1:100) with TSA-TYR-520 (Green), and E-Cadherin Rabbit mAb (A20798, 1:500) with TSA-TYR-570 (Red). DAPI (Blue) was used for nuclear staining. Prior to multiplex IHC staining, high-pressure antigen retrieval was performed using 0.01M citrate buffer at pH 6.0. The analysis was

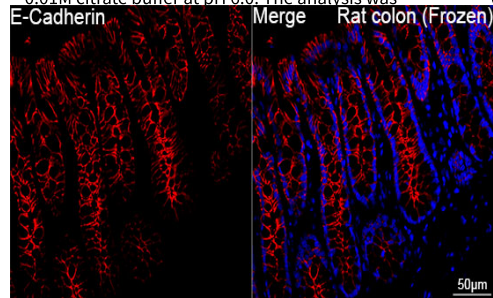

Confocal imaging of frozen sections of Rat colon tissue using E-Cadherin Rabbit mAb (A20798, dilution 1:200) followed by a further incubation with Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Microwave antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

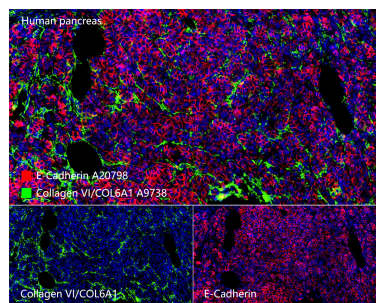

The multiplex IHC analysis on paraffin-embedded Human pancreas tissue using the following specific primary antibodies and tyramide signal amplification (TSA) reagents (RK05903) : Collagen VI/COL6A1 Rabbit mAb (A9738, 1:100) with TSA-TYR-520 (Green), and E-Cadherin Rabbit mAb (A20798, 1:500) with TSA-TYR-570 (Red). DAPI (Blue) was used for nuclear staining. Prior to multiplex IHC staining, high-pressure antigen retrieval was performed using 0.01M citrate buffer at pH 6.0. The analysis was completed using a 40x objective lens.

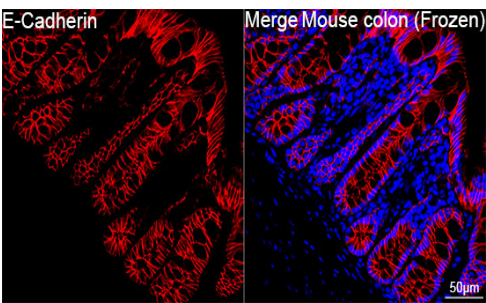

Confocal imaging of frozen sections of Mouse colon tissue using E-Cadherin Rabbit mAb (A20798, dilution 1:200) followed by a further incubation with Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Microwave antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

