

Cytokeratin 18 Rabbit mAb

Catalog NO.: A27298 **Recombinant**

Basic Information

Observed MW

48kDa

Calculated MW

48kDa

Category

Monoclonal Antibodies

Applications

WB, IF/ICC, IF-P, IHC-P, ELISA

Cross-Reactivity

Human, Mouse, Rat

Conjugation

Unconjugated

Clone number

ARC73118

Background

KRT18 encodes the type I intermediate filament chain keratin 18. Keratin 18, together with its filament partner keratin 8, are perhaps the most commonly found members of the intermediate filament gene family. They are expressed in single layer epithelial tissues of the body. Mutations in this gene have been linked to cryptogenic cirrhosis. Two transcript variants encoding the same protein have been found for this gene.

Recommended Dilutions

WB	1:12000 - 1:72000
IF/ICC	1:200 - 1:2000
IF-P	1:200 - 1:2000
IHC-P	1:6000 - 1:30000

ELISA

Recommended starting concentration is 1 µg/mL.
Please optimize the concentration based on your specific assay requirements.

Contact

Immunogen Information

Gene ID

3875

Swiss Prot

P05783

Immunogen

Recombinant protein (or fragment).This information is considered to be commercially sensitive.

Synonyms


K18; CK-18; CYK18

Product Information

☒		order@abclonal.com
☒		service@abclonal.com
☒		www.abclonal.com

Source	Isotype	Purification
Rabbit	IgG	Affinity purification
Storage		
Store at -20°C. Avoid freeze / thaw cycles.		
Buffer: PBS with 0.09% Sodium azide,0.05% BSA,50% glycerol,pH7.3.		

Validation Data

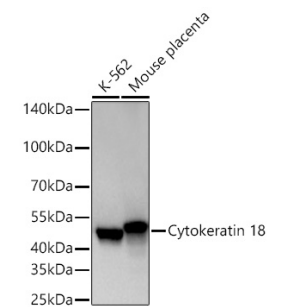

Western blot analysis of various lysates using Cytokeratin 18 Rabbit mAb (A27298) at 1:12000 dilution incubated overnight at 4°C.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25 µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 10s.

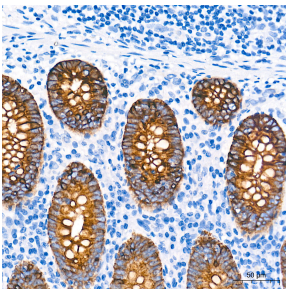

Immunohistochemistry analysis of paraffin-embedded Human appendix tissue using Cytokeratin 18 Rabbit mAb (A27298) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

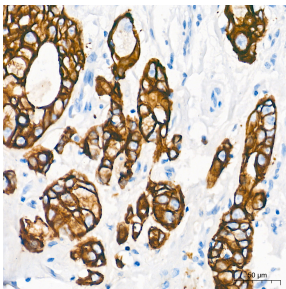


Immunohistochemistry analysis of paraffin-embedded Human breast cancer tissue using Cytokeratin 18 Rabbit mAb (A27298) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

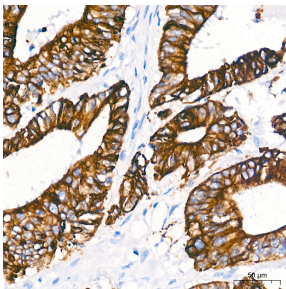


Immunohistochemistry analysis of paraffin-embedded Human colon carcinoma tissue using Cytokeratin 18 Rabbit mAb (A27298) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.



Immunohistochemistry analysis of paraffin-embedded Human fallopian tube tissue using Cytokeratin 18 Rabbit mAb (A27298) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

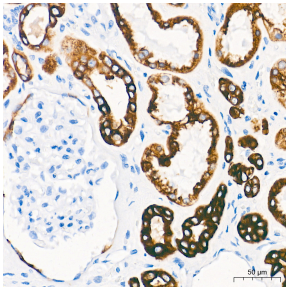

Immunohistochemistry analysis of paraffin-embedded Human kidney tissue using Cytokeratin 18 Rabbit mAb (A27298) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

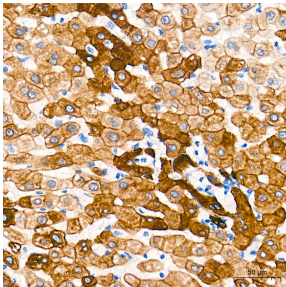


Immunohistochemistry analysis of paraffin-embedded Human liver tissue using Cytokeratin 18 Rabbit mAb (A27298) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.


High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.


Validation Data

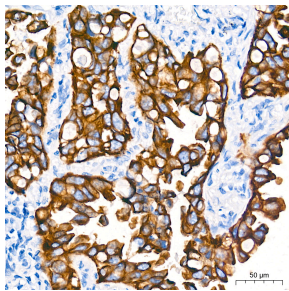


Immunohistochemistry analysis of paraffin-embedded Human lung adenocarcinoma tissue using Cytokeratin 18 Rabbit mAb (A27298) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

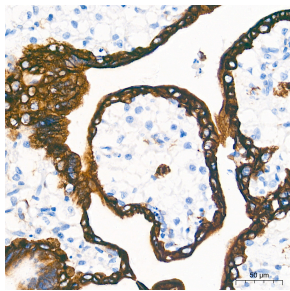


Immunohistochemistry analysis of paraffin-embedded Human placenta tissue using Cytokeratin 18 Rabbit mAb (A27298) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

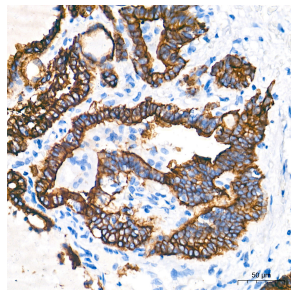


Immunohistochemistry analysis of paraffin-embedded Human thyroid cancer tissue using Cytokeratin 18 Rabbit mAb (A27298) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

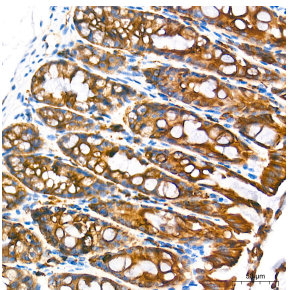


Immunohistochemistry analysis of paraffin-embedded Mouse colon tissue using Cytokeratin 18 Rabbit mAb (A27298) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

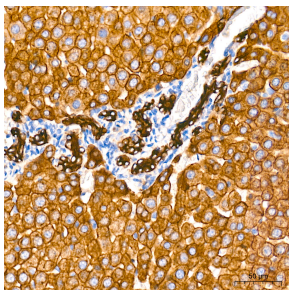


Immunohistochemistry analysis of paraffin-embedded Rat liver tissue using Cytokeratin 18 Rabbit mAb (A27298) at a dilution of 1:6000 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

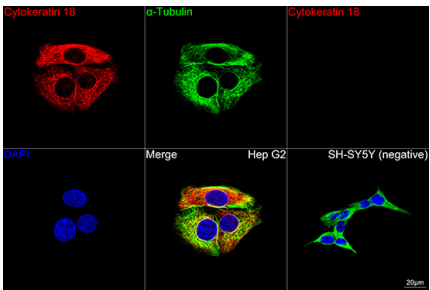


Confocal imaging of Hep G2 (positive) and SH-SY5Y (negative) cells using Cytokeratin 18 Rabbit mAb (A27298, dilution 1:800) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with α -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.

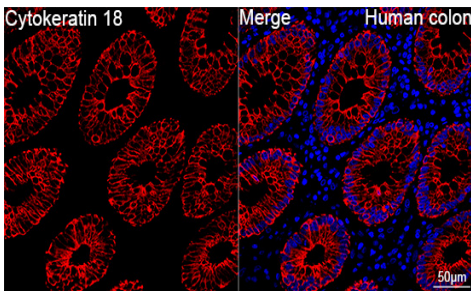

Confocal imaging of paraffin-embedded Human colon tissue using Cytokeratin 18 Rabbit mAb (A27298, dilution 1:800) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI

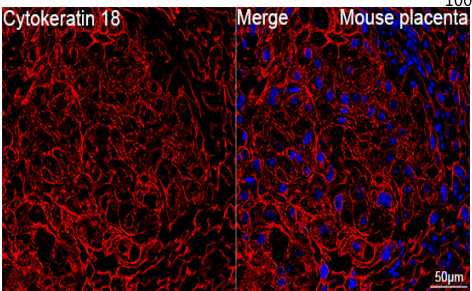

Confocal imaging of paraffin-embedded Mouse placenta tissue using Cytokeratin 18 Rabbit mAb (A27298, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI

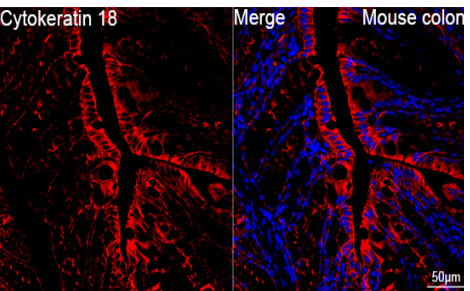


Confocal imaging of paraffin-embedded Mouse colon tissue using Cytokeratin 18 Rabbit mAb (A27298, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI

was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

Validation Data

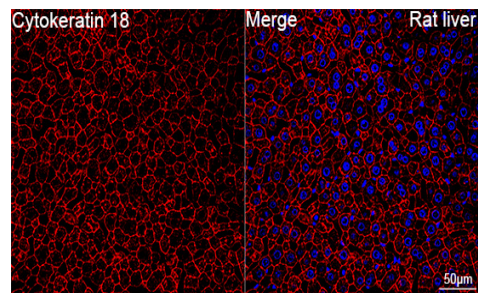

Confocal imaging of paraffin-embedded Rat liver tissue using Cytokeratin 18 Rabbit mAb (A27298, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

