

ESD Rabbit pAb

Catalog NO.: A6228

Basic Information

Observed MW

36kDa

Calculated MW

31kDa

Category

Polyclonal Antibodies

Applications

WB, IHC-P, ELISA

Cross-Reactivity

Human, Mouse

Conjugation

Unconjugated

Background

This gene encodes a serine hydrolase that belongs to the esterase D family. The encoded enzyme is active toward numerous substrates including O-acetylated sialic acids, and it may be involved in the recycling of sialic acids. This gene is used as a genetic marker for retinoblastoma and Wilson's disease.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:50 - 1:200

ELISA

Recommended starting concentration is 1 µg/mL.
Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID

2098

Swiss Prot

P10768

Immunogen

Recombinant protein (or fragment).This information is considered to be commercially sensitive.

Synonyms

FGH; ESD

Contact

Product Information

☒		order@abclonal.com
☒		service@abclonal.com
☒		www.abclonal.com

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage


Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.


Validation Data

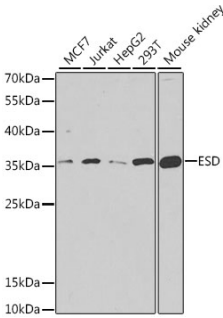

Western blot analysis of various lysates using ESD Rabbit pAb (A6228) at 1:1000 dilution.
 Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.
 Lysates/proteins: 25µg per lane.
 Blocking buffer: 3% nonfat dry milk in TBST.
 Detection: ECL Basic Kit (RM00020).
 Exposure time: 90s.

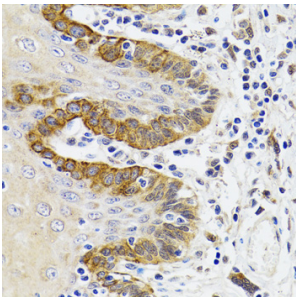

Immunohistochemistry analysis of paraffin-embedded Human trachea using ESD Rabbit pAb (A6228) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

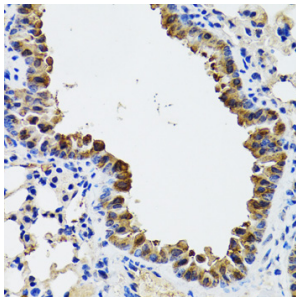

Immunohistochemistry analysis of paraffin-embedded Mouse lung using ESD Rabbit pAb (A6228) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

